

PRODUCT SPECIFICATIONS


Doctor C Dental Stool


www.australiandental.com.au/doctor-c-dental-stool




Tel: +61 3 9391 3552 Email: sales@australiandental.com.au Fax: +61 3 9923 6516 Factory address: 13/17 Chambers Road, Altona North, VIC 3025, Australia Postal address: PO Box 200, Altona North, VIC 3025, Australia

Doctor C Dental Stool

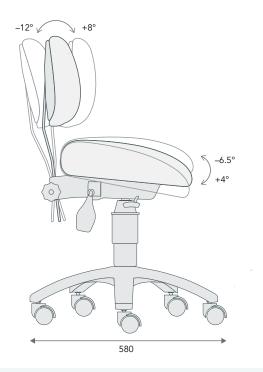

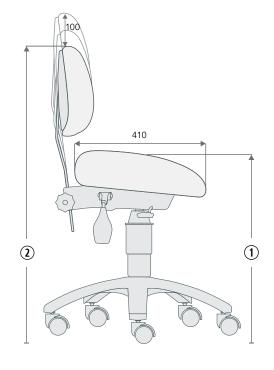

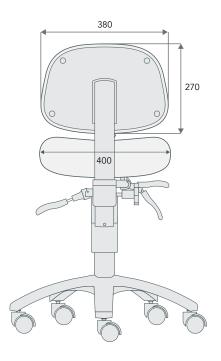

PRODUCT OPTIONS

The Doctor C Stool is fitted with a 140 stroke gas lift as standard. If you require a higher or lower stool, alternative gas lifts are available (see table below).

	Seat height $^{(1)}$		Back height ^②	
Gas lift size	MIN	МАХ	MIN	МАХ
100	470	545	710	795
140	490	610	735	860
160	545	680	790	930
200	550	740	800	990

Note:

- Heights are estimates only (tolerance 10mm)
- Vertical movement shown in diagram is in addition to height shown in table

- With alloy base fitted, all heights are reduced by 20mm

COMPONENTS

Dark grey Black White

Upholstery

The seat is manufactured with long-lasting injection moulded foam. Choose from over 200 vinyl colours or supply your own. For a full list of colours, see our website. All our vinyls are antibacterial and flame retardant.

Mechanism

Fully ergonomic 3 lever mechanism with independent seat and back adjustment.

Base material

Nylon (standard) or alloy (alloy diameter 600mm).

Friction castors

Fitted with polyurethane soft tire friction castors limiting the stool's movement when pressure is released (e.g. when the user stands up). Free wheel castors can be fitted upon request (black only).

CAPACITY

120kg (evenly spread).

SHIPPING INFO

Stool weight: 11.5kg Packaged weight: 13.5kg Packaged dimensions: 590mm × 590mm × 510mm

OPTIONAL ACCESSORIES

The height of the stool is normally controlled by the right hand lever. A foot activation ring can be fitted to all doctor stools (currently grey only).

HOW TO ORDER

Dental companies or government agencies, please contact us directly. Dentists please see our 'store locator' page for a list of our Australian and international distributors.