

**OptiDisc™** | High-performing polishing and contouring disc is now even easier to choose!



**FEATURES & BENEFITS**

**OptiDisc™**

- The translucent material affords great visibility of the working area, for extra precise finishing.
- Made of polyester coated with aluminium oxide particles.
- Patented mandrel coating for avoiding scratching or coloring of the teeth.
- Long-Lasting performance.
- 4-grit sizes for contouring, finishing, polishing and high gloss-polishing.
- The colour coded system allows you to easily identify the appropriate grit level based on your procedural needs.


**OptiDisc™ Mandrels**

- Patented mandrel coating for avoiding scratching or coloring of the teeth.
- Optimal torque transmission to disc, no sliding.
- Mandrels are autoclavable.
- 2 Mandrel sizes & lengths.
- The shorter mandrel facilitates an easy access to posterior and other restricted areas.

**Easy to recognise colour-coding**

Colour coded stages of abrasiveness means easy recognition of grit size and clinical step.


<b>CONTOUR</b>	<b>FINISH</b>	<b>POLISH</b>	<b>HIGH GLOSS POLISH</b>
<i>Extra-Coarse 80 µm</i>	<i>Coarse/Medium 40 µm</i>	<i>Fine 20 µm</i>	<i>Extra-Fine 10 µm</i>

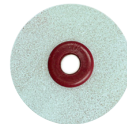

Active Side



Non Active Side

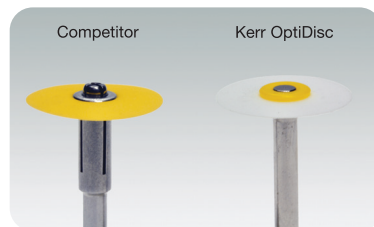

**Kerr Corporation**  
1717 W. Collins Ave.  
Orange, CA 92867 USA  
1-800-KERR-123  
kerrdental.com

**Kerr**

## Order Information

### OptiDisc

#### Assorted Kits

4200 **OptiDisc General Assorted Kit 9.6mm, 12.6mm (240 pcs)**

Contains:

∅ 9.6mm: 30/Extra-Coarse, 30/Coarse, 30/Fine, 30/Extra-Fine  
 ∅ 12.6mm: 30/Extra-Coarse, 30/Coarse, 30/Fine, 30/Extra-Fine  
 3 standard mandrels, 2 short mandrels, 1 OptiShine

4190 **OptiDisc Assorted kit 15.9mm (80 pcs)**

Contains:

∅ 15.9mm: 20/Extra-Coarse, 20/Coarse, 20/Fine, 20/Extra-Fine

4188 **OptiDisc Mini Kit 12.6mm (120 pcs)**

Contains:

∅ 12.6mm: 30/Extra-Coarse, 30/Coarse, 30/Fine, 30/Extra-Fine  
 1 standard mandrel, 1 short mandrel, 1 OptiShine

195 **OptiDisc standard mandrel (5 pcs)**

196 **OptiDisc short mandrel (5 pcs)**

#### Refill packs

##### ∅ 9.6mm

4197 **OptiDisc Extra-Coarse 9.6mm (80 pcs)**

4181 **OptiDisc Coarse/Medium 9.6mm (100 pcs)**

4182 **OptiDisc Fine 9.6mm (100 pcs)**

4183 **OptiDisc Extra-Fine 9.6mm (100 pcs)**

##### ∅ 12.6mm

4198 **OptiDisc Extra-Coarse 12.6mm (80 pcs)**

4184 **OptiDisc Coarse/Medium 12.6mm (100 pcs)**

4185 **OptiDisc Fine 12.6mm (100 pcs)**

4186 **OptiDisc Extra-Fine 12.6mm (100 pcs)**

##### ∅ 15.9mm

4199 **OptiDisc Extra-Coarse 15.9mm (80 pcs)**

4191 **OptiDisc Coarse/Medium 15.9mm (100 pcs)**

4192 **OptiDisc Fine 15.9mm (100 pcs)**

4193 **OptiDisc Extra-Fine 15.9mm (100 pcs)**

## Other related products:

### OptiStrip™

#### Finishing & Polishing Strip



#### FEATURES & BENEFITS

- Used for finishing and polishing of inter-dental surfaces.
- Ideal for polishing direct restorations including composites, glass ionomers, amalgams, semiprecious and precious metals.
- The abrasive particles have an optimized shape to reduce wear and keep sharpness for lasting efficiency.
- Offers added value in providing a faster finishing and polishing procedure, in only 2 steps.

#### Order Information

1990 **OptiStrip Assorted Kit**

Contains: 4 x 50 OptiStrip of each type

#### Refill Packs of 100

1991 **OptiStrip Finishing Strip, 2 mm**

1992 **OptiStrip Polishing Strip, 2 mm**

1993 **OptiStrip Finishing Strip, 4 mm**

1994 **OptiStrip Polishing Strip, 4 mm**

### Opti1Step™ Polisher

#### Composite One-Step High Gloss Polisher



#### FEATURES & BENEFITS

- Multiple use polishing system.
- Indicated for Composite.
- Increases Surface Gloss.
- Reduces surface roughness and bacterial retention.
- Cost-efficient compared to traditional procedures.
- Different shapes.
- Autoclavable at 134°C, reusable.
- Easy and fast technique: no paste needed.
- Gold-plated mandrel and unique matrix for greater durability.

#### Order Information

8000 **Opti1Step Assortment Kit**

Contains: 12 Composite Opti1Step Polisher: 4 Flames, 4 Cups, 4 Discs

8004 **Opti1Step Assorted Minikit**

Contains: 6 Composite Opti1Step Polisher: 2 Flames, 2 Cups, 2 Discs

8001 **Opti1Step** 

Contains: 12 Composite Polisher

8002 **Opti1Step** 

Contains: 12 Composite Polisher

8003 **Opti1Step** 

Contains: 12 Composite Polisher

### OptiShine™

### Occlubrush™

#### Polishing Brushes



#### FEATURES & BENEFITS

##### OptiShine™

- The first concave-shaped polishing brush.
- Autoclavable for multiple units, at 134°C, at least 3 min. No effects on polishing performance.
- Increases surface gloss.
- Reduces roughness and bacterial retention.

##### Occlubrush™

- Non-destructive polishing of occlusal and concave surfaces.
- Past free.
- Faster and more efficient than polishing paste systems.
- The bristles' special fiber allows a high gloss shine.
- 3 shapes to reach all areas and angles.

#### Order Information

**OptiShine™**

pack of 3    pack of 10  
2513    2514

#### Occlubrush™

2520 **Occlubrush Assorted Kit**

Contents: 5 Occlubrush regular cup, 5 Occlubrush small cup, 5 Occlubrush point

#### Refill Packs

pack of 3    pack of 10

**Occlubrush™ Regular Cup** 2503    2510

**Occlubrush™ Small Cup** 2504    2511

**Occlubrush™ Point** 2505    2512